

Perfect Ergonomic System

FlexFit®




Flexible Proximal End

- Requires less insertion force than competitors' to make a secure seal, reducing the risk of repetitive stress injury
- Tips don't fall off multichannel pipettes
- Completes the ergonomic system on LTS pipettes

StarStop®

Positive Stop Feature

- Gives consistently low ejection forces on LTS pipettes
- Aligns the tips for even liquid levels on multichannels

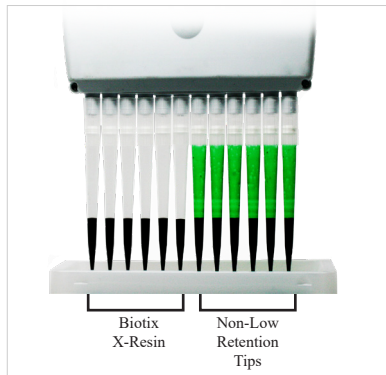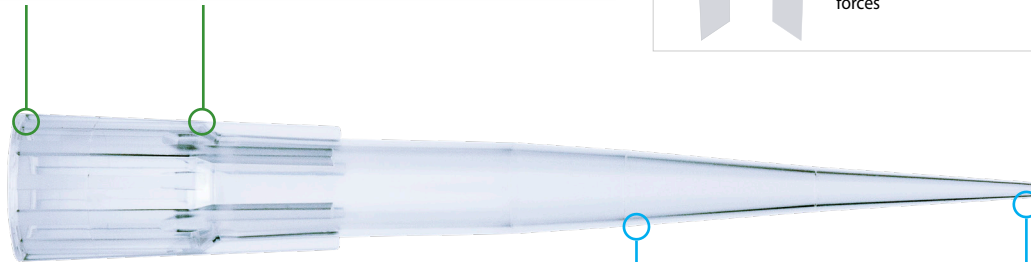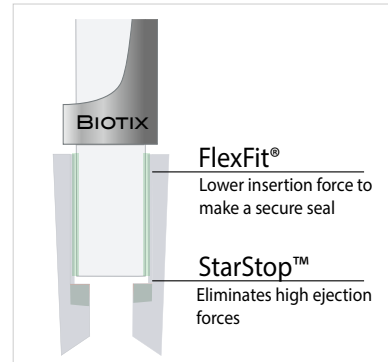

Clean Delivery

X-Resin®

Low Retention Surface

- Maximizes sample recovery
- Improves accuracy across your experiments

Blade®

Minimizes Surface Tension

- Eliminates hanging droplet formation
- Saves time by reducing the need for tip touch off
- Increases accuracy and precision

Biotix® xTIP4 Pipette Tips, Low Retention

Volume	20 µL		200 µL	250 µL		1000 µL	
Tips/Pack	960		960	960		768	
Sterilized	No	Yes	Yes	No	Yes	No	Yes
Racked, Filtered		63300000	63300001				63300003
Racked	63300020	63300005		63300021	63300006	63300023	63300008
CleanPak Reload	63300025	63300010		63300026	63300011	63300028	63300013
Tip Eject Reload		63300015			63300016		63300018

xTIP4™

Rainin LTS Compatible Pipette Tips

Optimus Rack	
Structural Support	<ul style="list-style-type: none"> Rigid lid and base design for repeatable autoclaving Stable base profile inhibits tipping of the rack when pipetting the corners
Locking Clasp	<ul style="list-style-type: none"> Secures tip during transport Clicks to ensure rack is securely closed
Snap Card Tabs	<ul style="list-style-type: none"> Simply squeeze and pull up for easy one handed snap card removal with reloads
Double Ball Hinge	<ul style="list-style-type: none"> Easy to remove and reattach Opens past 180° for optimal multichannel pipetting
Modified Hole Design	<ul style="list-style-type: none"> Prevents rack lift up when removing tips at an angle



Tip Eject Reload		
Convenience <ul style="list-style-type: none"> Clear casing and colored snap cards for easy tip volume identification Snap cards click into place Easily operated with one hand without contact with tips or snap card 	Sustainability <ul style="list-style-type: none"> 100% recyclable packaging Empty racks can be autoclaved for multiple uses 	Contamination Prevention <ul style="list-style-type: none"> Tips are shielded from the environment and user handling Rigid plastic case provides a cleaner reload than cardboard alternatives All Tip Eject reloads come pre-sterilized and ready to use

CleanPak Reload		
Convenience <ul style="list-style-type: none"> CleanPak uses 1/5 the footprint as a pack of racked tips Inserts are easy to transfer to your empty rack, one at a time or all at once Clear bag and colored snapcards make tip volume identification easy 	Sustainability <ul style="list-style-type: none"> 100% recyclable packaging Uses ~50% less plastic than racked tips 	Contamination Prevention <ul style="list-style-type: none"> Resealable bag with sterilized and non-sterilized options Arrives hermetically sealed Contains no contaminating cardboard or paper



Biotix Inc
9880 Mesa Rim Road
San Diego, CA 92121
United States

Customer Service: +1 858-875-7696 Monday-Friday 6AM - 5PM (PST)

To Contact a sales rep visit: www.biotix.com/contact

Biotix
Fluid Innovation