

## Perfect Ergonomic System

### FlexFit®




#### Flexible Proximal End

- Requires less insertion force than competitors' to make a secure seal, reducing the risk of repetitive stress injury
- Tips don't fall off multichannel pipettes
- Completes the ergonomic system on LTS pipettes

### StarStop™

#### Positive Stop Feature

- Gives consistently low ejection forces on LTS pipettes
- Aligns the tips for even liquid levels on multichannels

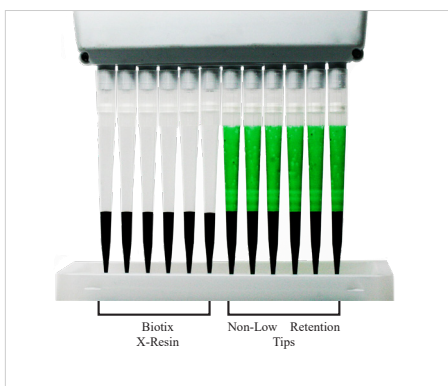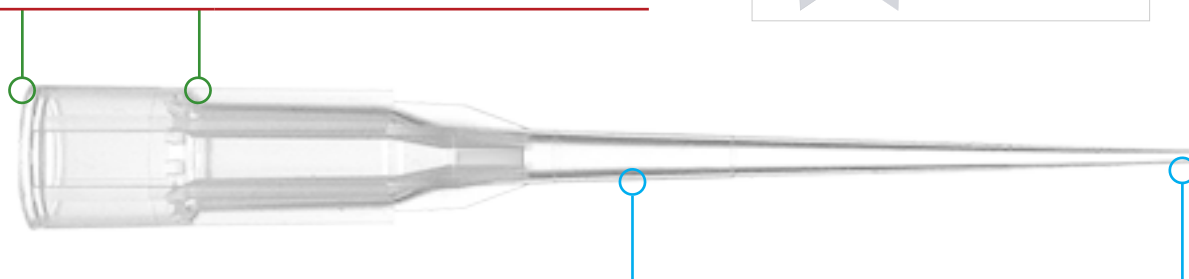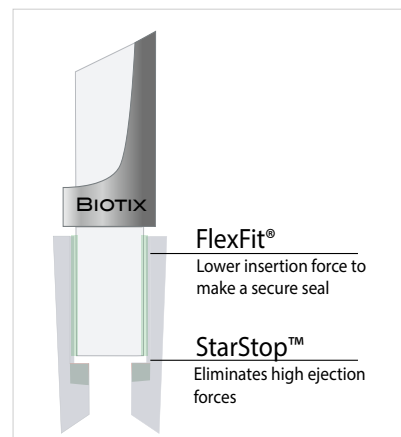

## Clean Delivery

### X-Resin®

#### Low Retention Surface

- Maximizes sample recovery
- Improves accuracy across your experiments

### Blade®

#### Minimizes Surface Tension

- Eliminates hanging droplet formation
- Saves time by reducing the need for tip touch off
- Increases accuracy and precision

Volume $\mu$ L	Pre-Sterilized	Racked, Filtered	Racked	CleanPak Reload	Tip Eject Reload	Packaging
20	Yes	63300931	63300005	63300025	63300015	960 tips/pack
	No		63300020			
200	Yes	63300001				960 tips/pack
	No					
250	Yes		63300006	63300026	63300016	960 tips/pack
	No		63300021			
300	Yes	63300002	63300007	63300027		960 tips/pack
	No		63300022			
1000	Yes	63300003	63300008	63300028	63300018	768 tips/pack
	No		63300023			
1200	Yes	63300004	63300009	63300029	63300019	768 tips/pack
	No		63300024			

# A More Sustainable Rack

xTIP

The Biotix Rack provides top functionality and strength while using 1/3 less plastic than leading competitors



A more sustainable design with all the desired features:

- More rigid lid and base design for repeatable autoclave
- Locking clasp clicks to ensure rack is securely closed during transport
- Snap cards are easily removed with reloads by one hand
- Modified hole design prevents rack lift up during tip removal at any angle
- Lid opens past 180° and is easily removeable for optimal multichannel pipetting
- Compatible with both CleanPak and Tip Eject Reloads

## Reload Solutions

Biotix xTIP4 has two reload styles to reduce plastic waste in your lab. Simply reuse your empty racks.

### xTIP Tip Eject™ Reload



Tip Eject is a self aligning reload solution. Simply press on the top of the stack and the empty rack will be refilled at the bottom.

- Reduces waste by 30% when compared to racked product
- Sterilized after packaging for added protection
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors

### xTIP CleanPak™ Reload



CleanPak reloads arrive hermetically sealed in a ziploc bag for superior tip purity. Available in both sterilized and non-sterile formats.

- Reduces waste by 50% when compared to racked product
- Color coded snapcards match the compatible pipette's plunger cap
- To refill an empty rack just reach into the bag and manually transfer a snapcard
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors