

# Biotix xTIP Pipette Tips

xPIPETTE™ and LTS™ Compatible

## Perfect Ergonomic System

### FlexFit®




#### Flexible Proximal End

- Requires less insertion force than competitors' to make a secure seal, reducing the risk of repetitive stress injury
- Tips don't fall off multichannel pipettes
- Completes the ergonomic system on LTS pipettes

### StarStop™

#### Positive Stop Feature

- Gives consistently low ejection forces on LTS pipettes
- Aligns the tips for even liquid levels on multichannels

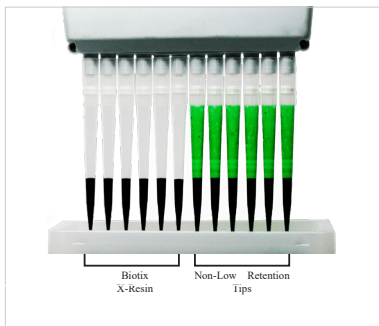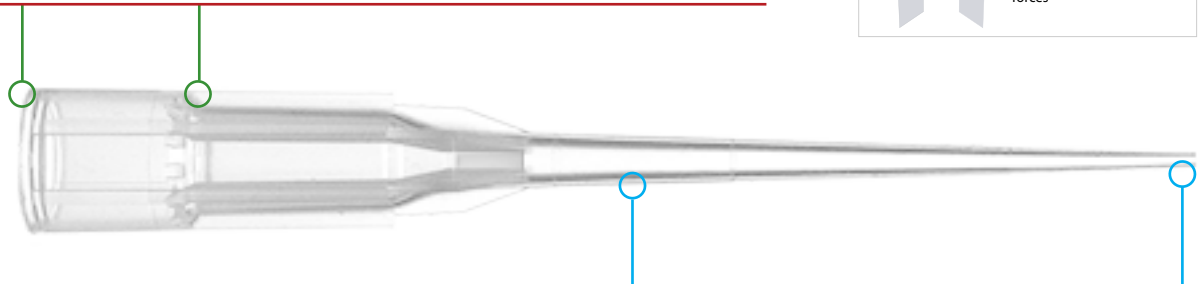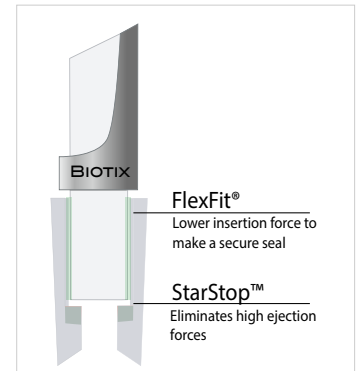

## Clean Delivery

### X-Resin®

#### Low Retention Surface

- Maximizes sample recovery
- Improves accuracy across your experiments

### Blade®

#### Minimizes Surface Tension

- Eliminates hanging droplet formation
- Saves time by reducing the need for tip touch off
- Increases accuracy and precision

Volume $\mu$ L	Pre-Sterilized	Racked, Filtered	Racked	CleanPak Reload	Tip Eject Reload	Packaging
20	Yes	76314-472	76266-150		76266-128	960 tips/pack
	No		76266-098	76266-118		
200	Yes	76266-142				960 tips/pack
	No					
250	Yes		76266-152		76266-130	960 tips/pack
	No		76266-100	76266-120		
300	Yes	76266-144	76266-092			960 tips/pack
	No		76266-102	76266-122		
1000	Yes	76266-146	76266-094		76266-134	768 tips/pack
	No		76266-104	76266-124		
1200	Yes	76266-148	76266-096		76266-136	768 tips/pack
	No		76266-106	76266-126		

# A More Sustainable Rack

The Biotix Rack provides top functionality and strength while using 1/3 less plastic than leading competitors



A more sustainable design with all the desired features:

- More rigid lid and base design for repeatable autoclave
- Locking clasp clicks to ensure rack is securely closed during transport
- Snap cards are easily removed with reloads by one hand
- Modified hole design prevents rack lift up during tip removal at any angle
- Lid opens past 180° and is easily removeable for optimal multichannel pipetting
- Compatible with both CleanPak and Tip Eject Reloads

## Reload Solutions

Biotix xTIP has two reload styles to reduce plastic waste in your lab. Simply reuse your empty racks.



### xTIP Tip Eject™ Reload

Tip Eject is a self aligning reload solution. Simply press on the top of the stack and the empty rack will be refilled at the bottom.

- Reduces waste by 30% when compared to racked product
- Sterilized after packaging for added protection
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors



### xTIP CleanPak™ Reload

CleanPak reloads arrive hermetically sealed in a ziploc bag for superior tip purity. Available in both sterilized and non-sterile formats.

- Reduces waste by 50% when compared to racked product
- Color coded snapcards match the compatible pipette's plunger cap
- To refill an empty rack just reach into the bag and manually transfer a snapcard
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors