

Perfect Ergonomic System

FlexFit®




Flexible Proximal End

- Requires less insertion force than competitors' to make a secure seal, reducing the risk of repetitive stress injury
- Tips don't fall off multichannel pipettes
- Completes the ergonomic system on LTS pipettes

StarStop™

Positive Stop Feature

- Gives consistently low ejection forces on LTS pipettes
- Aligns the tips for even liquid levels on multichannels

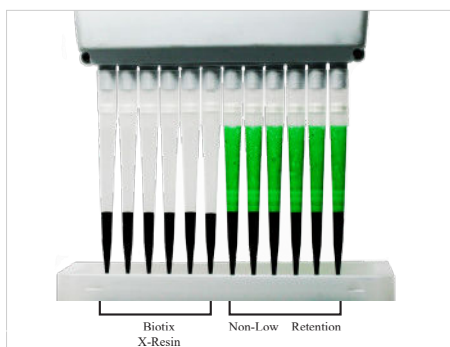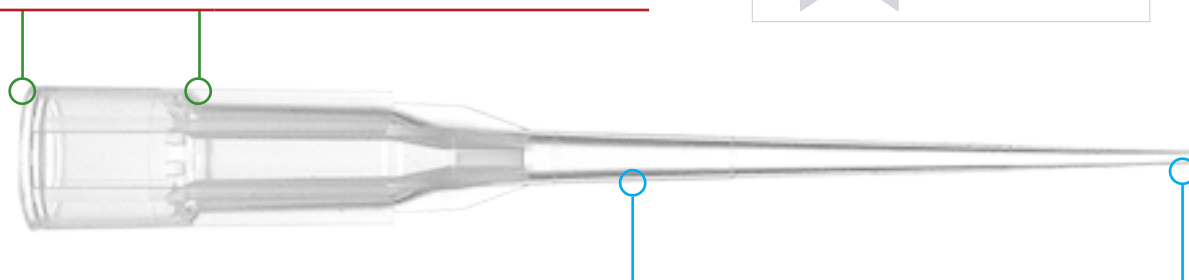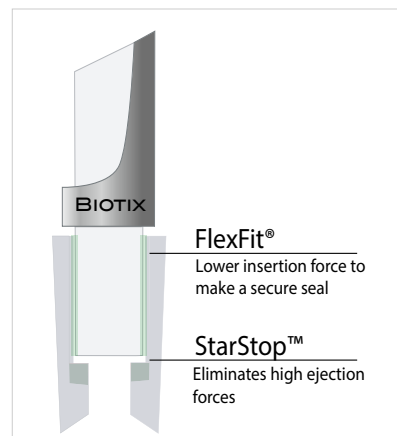

Clean Delivery

X-Resin®

Low Retention Surface

- Maximizes sample recovery
- Improves accuracy across your experiments

Blade®

Minimizes Surface Tension

- Eliminates hanging droplet formation
- Saves time by reducing the need for tip touch off
- Increases accuracy and precision

| Volume μ L | Sterilized | Racked, Filtered | Racked | CleanPak Reload | Tip Eject Reload | Packaging |
|----------------|------------|------------------|-----------|-----------------|------------------|---------------|
| 20 | Yes | 1156T51 | 1156T56 | 1156T71 | 1156T76 | 960 tips/pack |
| | No | | 1156T61 | | | |
| 200 | Yes | 1156T52 | | | | 960 tips/pack |
| | No | | | | | |
| 250 | Yes | | 1156T57 | 1156T72 | 1156T77 | 960 tips/pack |
| | No | | 1156T62 | | | |
| 300 | Yes | 23A00H173 | 23A00H174 | | | 960 tips/pack |
| | No | | 63300022 | | | |
| 1000 | Yes | 1156T54 | 1156T59 | 1156T74 | 1156T79 | 768 tips/pack |
| | No | | 1156T64 | | | |
| 1200 | Yes | 1156T55 | 1156T60 | 1156T75 | | 768 tips/pack |
| | No | | 1156T65 | | | |

A More Sustainable Rack

xTIP

The Biotix Rack provides top functionality and strength while using 1/3 less plastic than leading competitors



A more sustainable design with all the desired features:

- More rigid lid and base design for repeatable autoclave
- Locking clasp clicks to ensure rack is securely closed during transport
- Snap cards are easily removed with reloads by one hand
- Modified hole design prevents rack lift up during tip removal at any angle
- Lid opens past 180° and is easily removeable for optimal multichannel pipetting
- Compatible with both CleanPak and Tip Eject Reloads

Reload Solutions

Biotix xTIP4 has two reload styles to reduce plastic waste in your lab. Simply reuse your empty racks.

xTIP Tip Eject™ Reload



Tip Eject is a self aligning reload solution. Simply press on the top of the stack and the empty rack will be refilled at the bottom.

- Reduces waste by 30% when compared to racked product
- Sterilized after packaging for added protection
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors

xTIP CleanPak™ Reload



CleanPak reloads arrive hermetically sealed in a ziploc bag for superior tip purity. Available in both sterilized and non-sterile formats.

- Reduces waste by 50% when compared to racked product
- Color coded snapcards match the compatible pipette's plunger cap
- To refill an empty rack just reach into the bag and manually transfer a snapcard
- Rack boxes are autoclaveable after loading at 121° C
- Certified free of DNase, RNase, endotoxins, pyrogens, ATP, and PCR inhibitors